

BellaVitaTM

I N J E C T O R

The safer, easier and
faster drug delivery
system



For Better Results

Crow's Feet
Perioral Area
Acne Scars
Decolletage
Neck Lines
Forehead
Cellulite
Temple
Hands
Pores
Thigh



BellaVitaTM System

I N J E C T O R


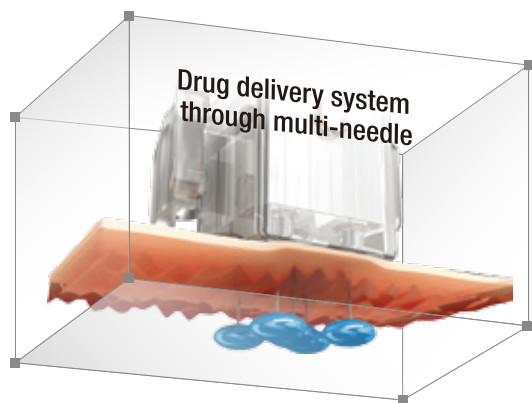
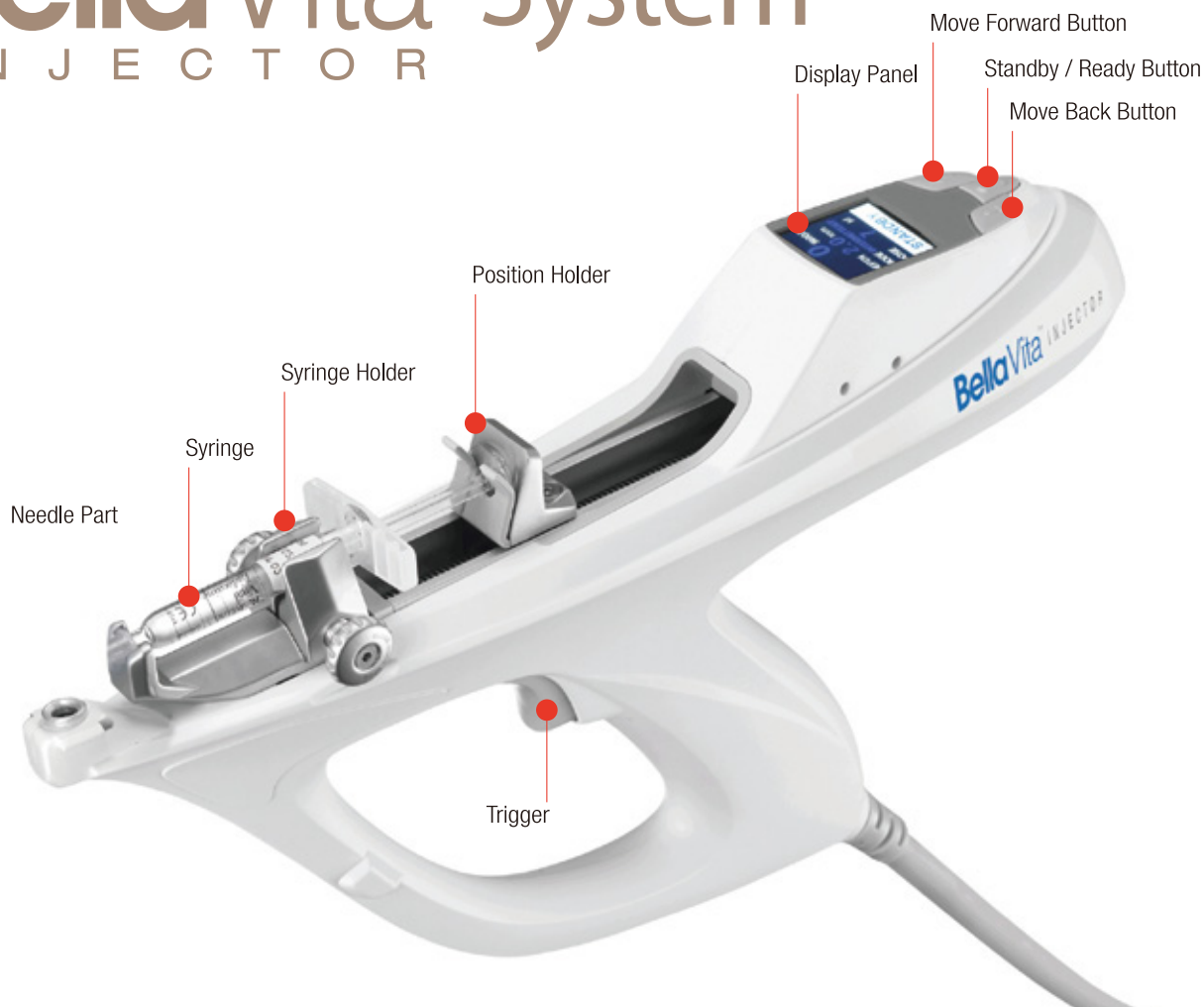


BellaVita InjectorTM is a specialized system for drug-delivery treatment.

The system is applied to the precise mechanical control, safety function and vacuum technology to obtain a stable treatment result.


There is 5-pin multi-needle to maximize the effectiveness of drug delivery in the target area. The system also offers advanced user convenient interface, users can easily operates main functions so that the skin rejuvenation treatment can be completed.

SKIN REJUVENATION

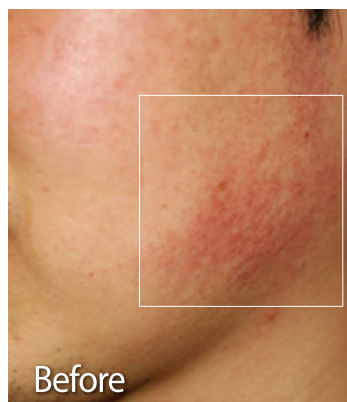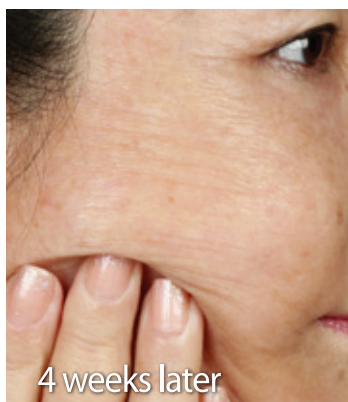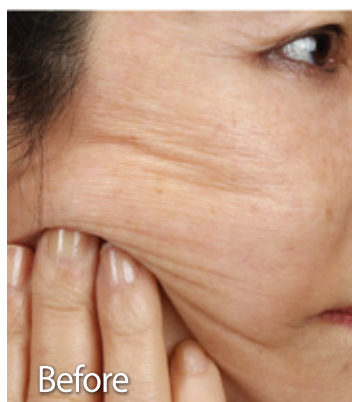

BellaVitaTM is your guide to more Professional treatment.

I N J E C T O R

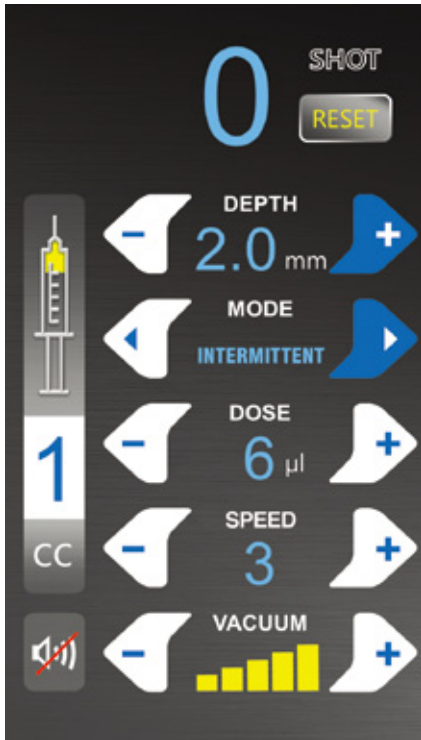

User Convenient Program

- Applied core functions
- Designed to be extremely easy to use

Minimizes Drug Loss

- Minimizes drug loss structure of multi-needle
- "SPEED" Function
minimize drug loss of drug viscosity and treatment dosage



Safe Treatment Available

- Provides stable mechanical control system
- Accurate needling and injection
- Quick and correct target area treatment available

Accurate Treatment Available

Vacuum function

- Guides stable treatment
- Reduces drug loss
- Minimizes bleeding
- Reducing side effect risk

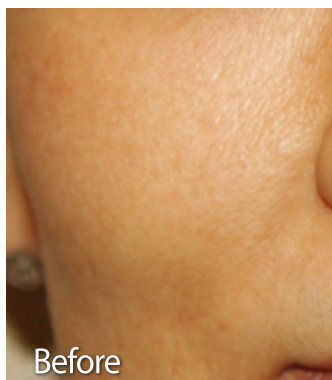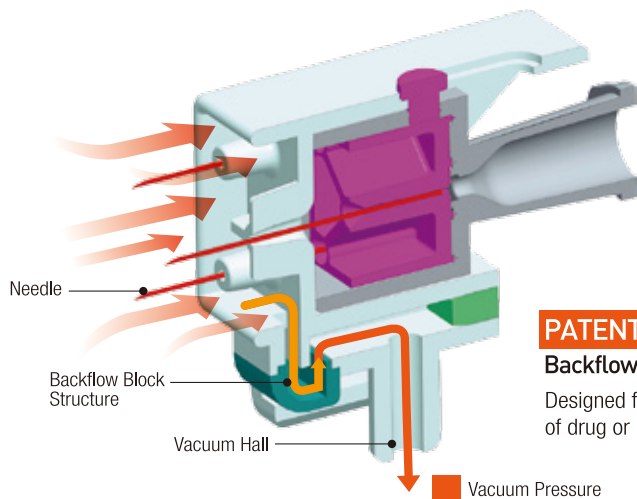

BellaVitaTM

I N J E C T O R




Installation

1. Turn on the power




Turn on the power switch.

2. Connecting multi-needle

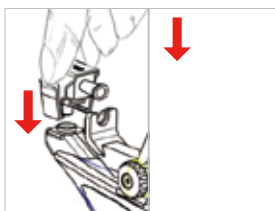

Match the multi-needle vacuum hall to the injector vacuum hall, and then press the upper part of the multi-needle until you hear a "click" sound.

3. Remove safe tape

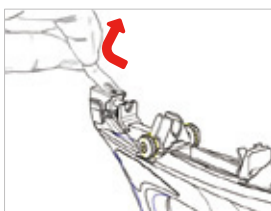


Remove the safe tape on the multi-needle.

6. Fasten the syringe holder




Fasten the syringe holder.

5. Connecting syringe and multi-needle

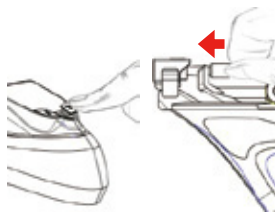


Connect the syringe to the multi-needle connect hall by pressing the forward button. After that push the syringe tightly by hand.

4. Syringe Installation

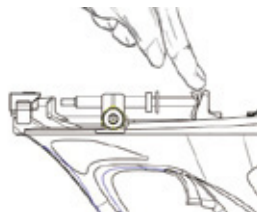

Put the syringe on the injector, then place the syringe piston correctly in the position holder.

Preparation

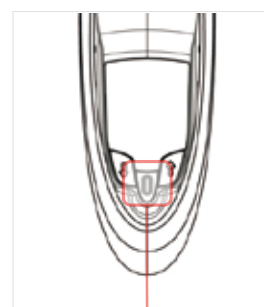

Standby/Ready Button

- Standby : Treatment standby mode
- Ready : Treatment ready mode



Move Back Button

- Move back the position holder.
- When pressed for more than two seconds, the position holder move back automatically, and press shortly once is to stop the motion.



Move Forward Button

- Move forward the position holder.

SPECIFICATION

Product Name	Bellavita Injector
Needle Penetration Depth	0~5 mm (interval 0.2mm)
Syringe Type	1cc, 2cc, 3cc, 5cc
Operation Mode	Intermittent, Continuous
Voltage & Frequency	DC 12V or DC 24V / 50~60Hz
Dimension (WxDxH mm)	Injector: 310 x 52 x 172 Body : 232 x 248 x 128
Weight	2.4 Kg (Body + Injector)

